

GC/MS Application Note



FOOD SAFETY



ENVIRONMENTAL

**Profiling of Gin Botanicals by
Headspace-SPME-ARROW-GCXGC-MSD/FID**



www.palsystem.com



Profiling of Gin Botanicals by Headspace-SPME-ARROW-GCXGC-MSD/FID

Andrew Ward, JSB® UK & Ireland 2018

Introduction

Gins are essentially neutral spirit flavoured with various botanicals by compounding or by distillation. There are hundreds of botanicals used to flavour gin, most commonly juniper with other seeds, herbs, roots and fruits. In distilled gins, the aromatic compounds (usually oils) are absorbed by the spirit in which the botanicals are steeped, or by the vapour as it passes through, so flavouring the distillate. The large numbers of compounds contributing to the aroma and flavour of gins coupled with the fact that flavour- and fragrance-key components can be present at very low concentrations makes chemical profiling challenging.

GCxGC-MS is the ideal analytical technique to accommodate the complexity of the samples due to its enhanced peak capacity and sensitivity. For quantitative purposes, and when using flow modulation, splitting to a FID is a logical choice. Sample preparation and sample introduction introduce further challenges due to the complex nature of the samples. Headspace-SPME Arrow is a fully automated sample preparation technique with the ability to extract targeted key compounds present at low concentrations whilst avoiding matrix compounds present at high concentrations. The combination of HSPME Arrow with GCxGC-MSD/FID is a very powerful tool in the comparative analysis of samples with complex matrices.

Samples

A set of twelve gins produced using the same base spirit was purchased from www.gininatin.co.uk. Each gin was described on its label with three of its major botanicals.

Instrumentation and Software

Analysis was fully automated using a CTC PAL RTC set up for the SPME Arrow technique, mounted on an Agilent 7890B GC equipped with an Agilent reverse-low modulator, FID and 5977B MSD with HES ion source. PAL Sample Control software

was used for the automated analysis. This software controls the Masshunter software which runs in the background. Zoex GCImage software was used for 2D data analysis.

Automation Workflow

The SPME Arrow technique is summarised in figure 1 and involves the following steps:

1. Incubation of sample in the (heated) agitator module
2. SPME Arrow Heatex Extraction in the heatex stirrer module
3. Desorption in the GC injector and triggering of the GC analysis
4. Conditioning of the Arrow in the SPME Arrow conditioning module

A new, triple-phase DVB/Carbon WR/PDMS Arrow was used for the analysis in order to extract as wide a range of target analytes as possible.

A video of the fully automated SPME Arrow technique can be found here:

<https://youtu.be/pEB7qT1TUww>

Figure 2 shows that this technique provides good resolution between compounds which are not easily resolved in 1D-GC. This is further illustrated in figure 3, which shows a 1D view taken from the same sample. All samples generated complex chromatographic profiles. Figure 4 shows that the samples can be clearly distinguished.

Two of the gin samples shared two of their three major botanicals - namely Orange Peel and Nutmeg, with one sample also listing Allspice and the other listing Pumpkin on the label. Whilst some similarities can be observed in Figure 4 between these samples in terms of the compounds present, the relative amounts of the majority of these compounds are clearly different. Furthermore, the 3D plot for the sample containing Allspice is dominated by the Estragole peak (see Figure 5).

PAL SPME Arrow: The Better SPME

PAL SYSTEM
Ingenious sample handling

Fully automated process → process safety and productivity

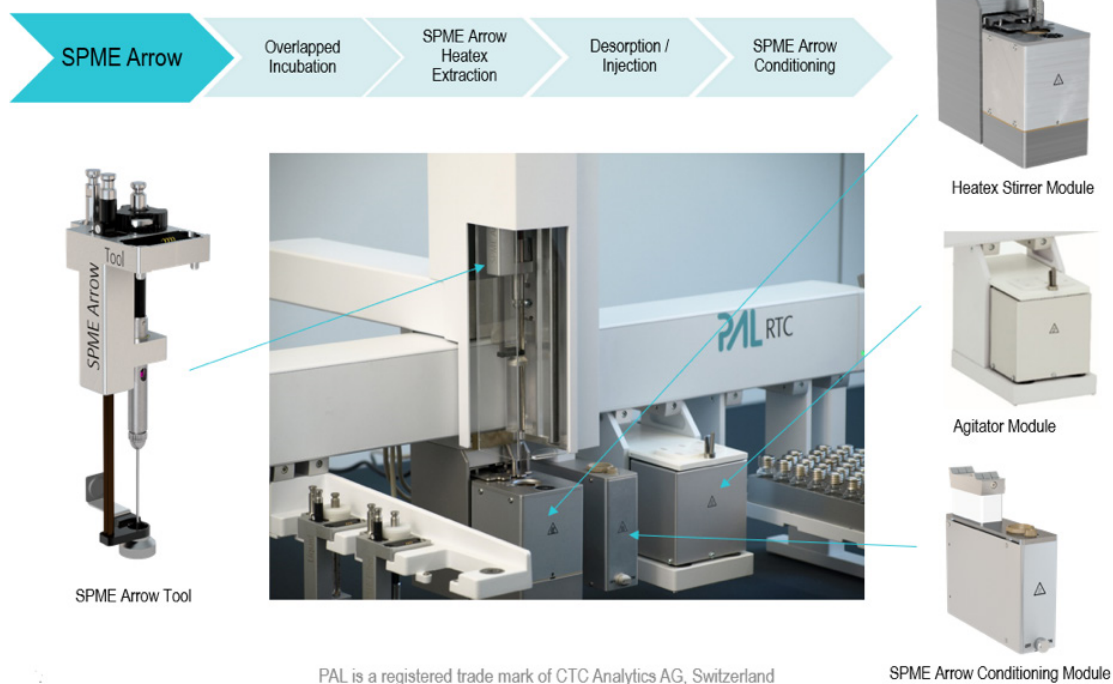


Fig. 1: The SPME Arrow technique on the CTC PAL RTC instrument.

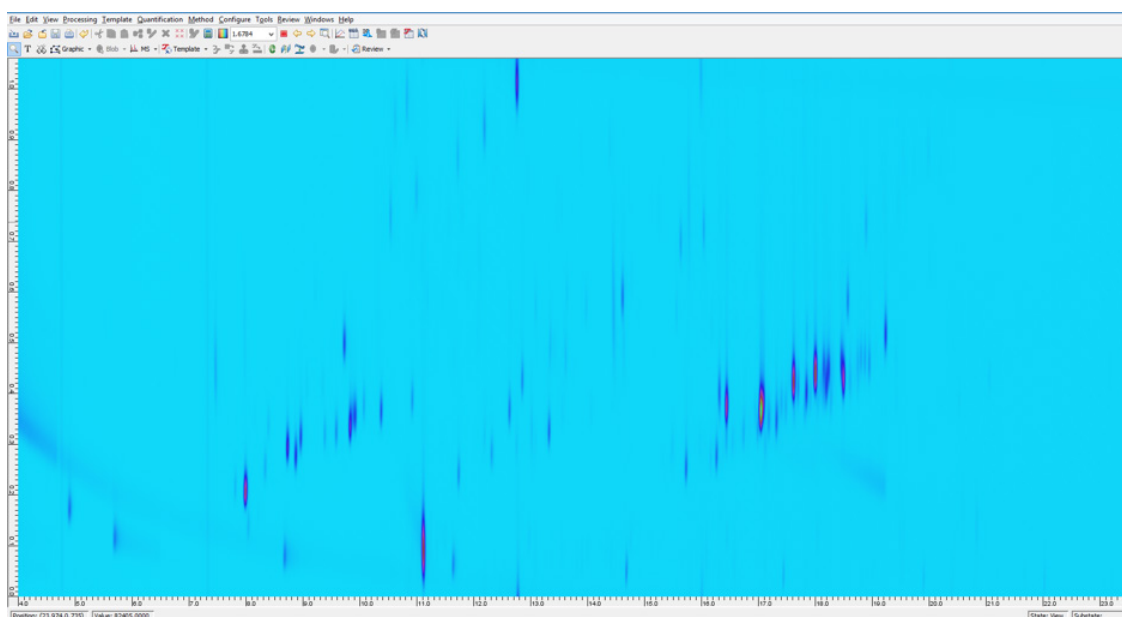


Fig. 2: 2D plot of sample "Orange Peel, Pumpkin, Nutmeg" (FID).



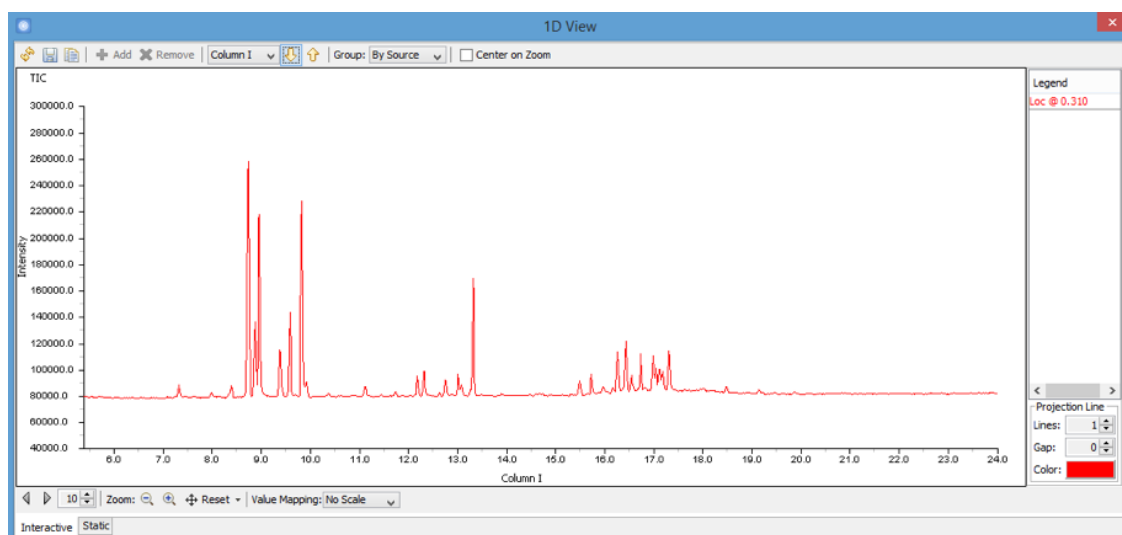



Fig. 3: Example 1D pot showing unresolved peaks in sample "Orange Peel, Pumpkin, Nutmeg".

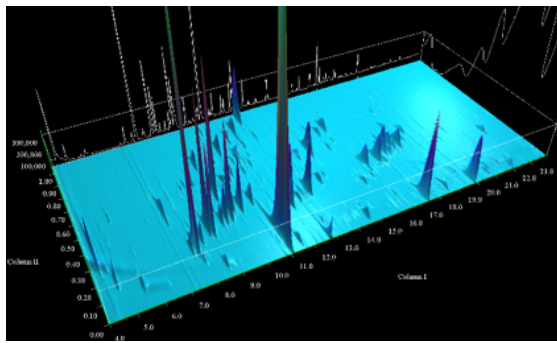


Sample "Angelica Root, Lemon Peel, Ginger"

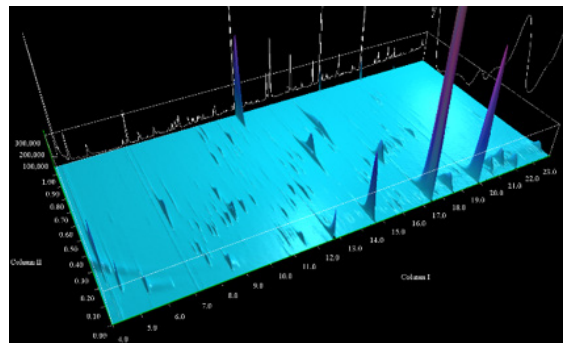


Sample "Black Pepper, Cassia Bark, Cloves"

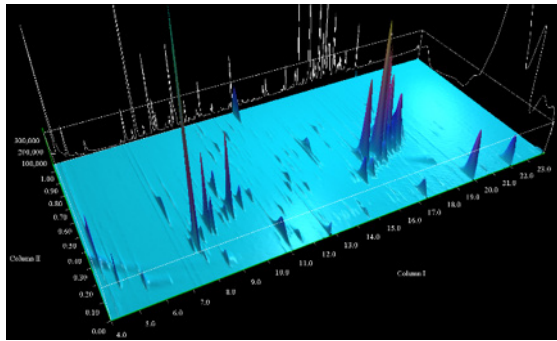


Sample "Cinnamon, Rhubarb, Ginger"

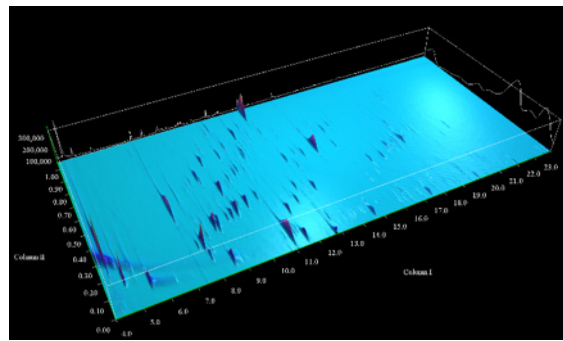


Sample "Cardamom, Samphire, Angelica"

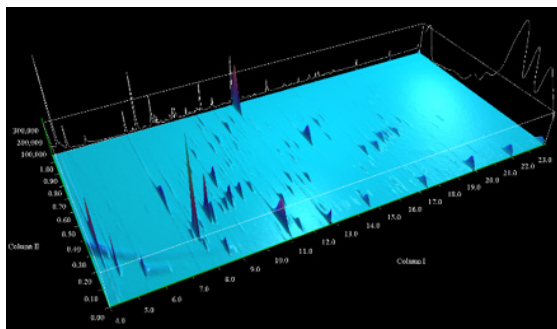


Sample "Elderflower, Orris Root, Apricot"

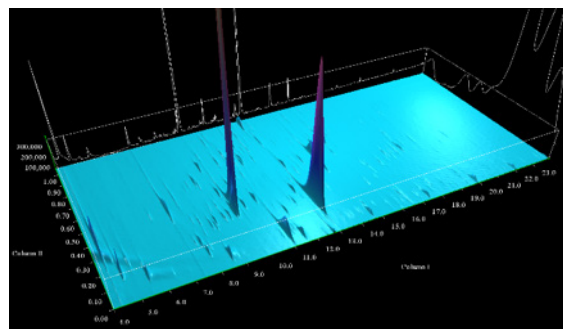


Sample "Fennel, Lemon, Thyme"

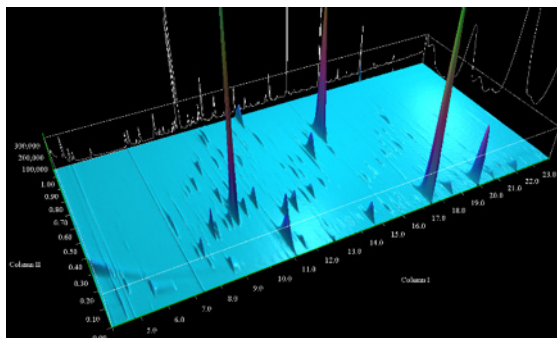


Sample "Lemon Peel, Coriander, Cardamom"

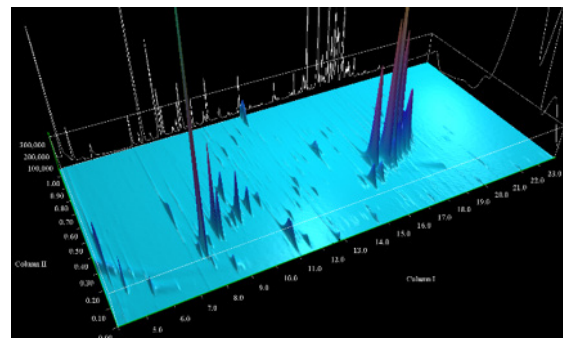

Sample "Maris Otter Malt, Coriander Seed, Fuggles Hops"

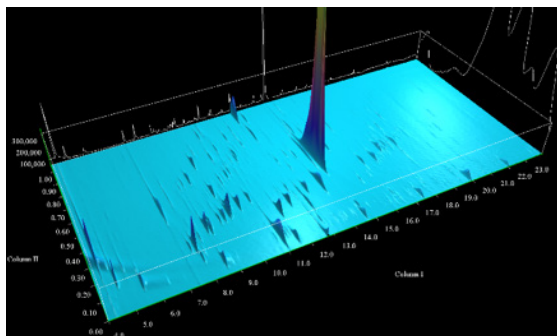


Sample "Orange Peel, Allspice, Nutmeg"

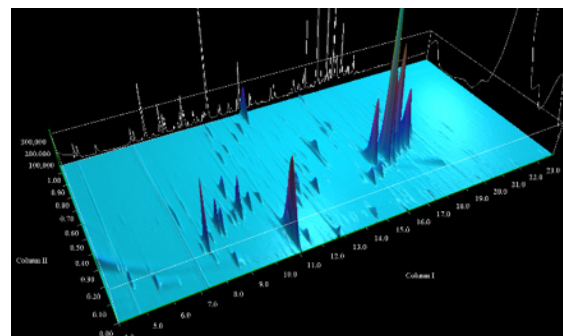


Sample "Orange Peel, Allspice, Nutmeg"

Fig. 4: 3D-plots of the 12 gin samples analyzed (FID).

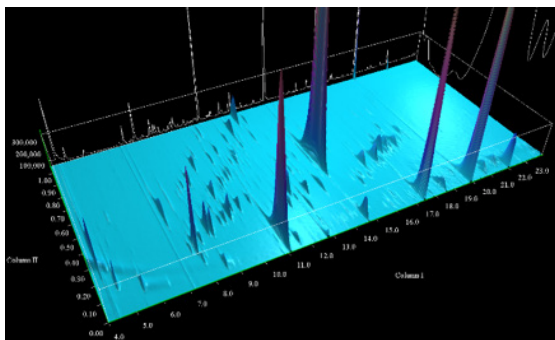

Sample "Packham Pear, Orange Peel, Cinnamon"

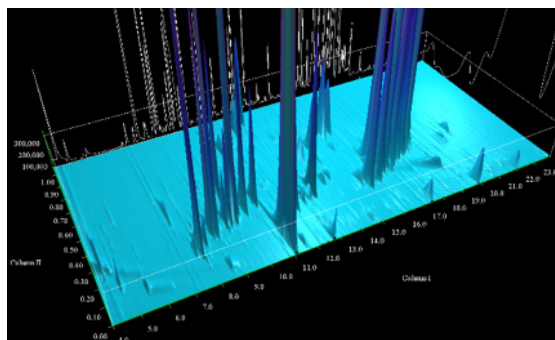

Sample "Pomegranate, Raspberry, Cardamom"

Fig. 4: 3D-plots of the 12 gin samples analysed (FID). (continued)

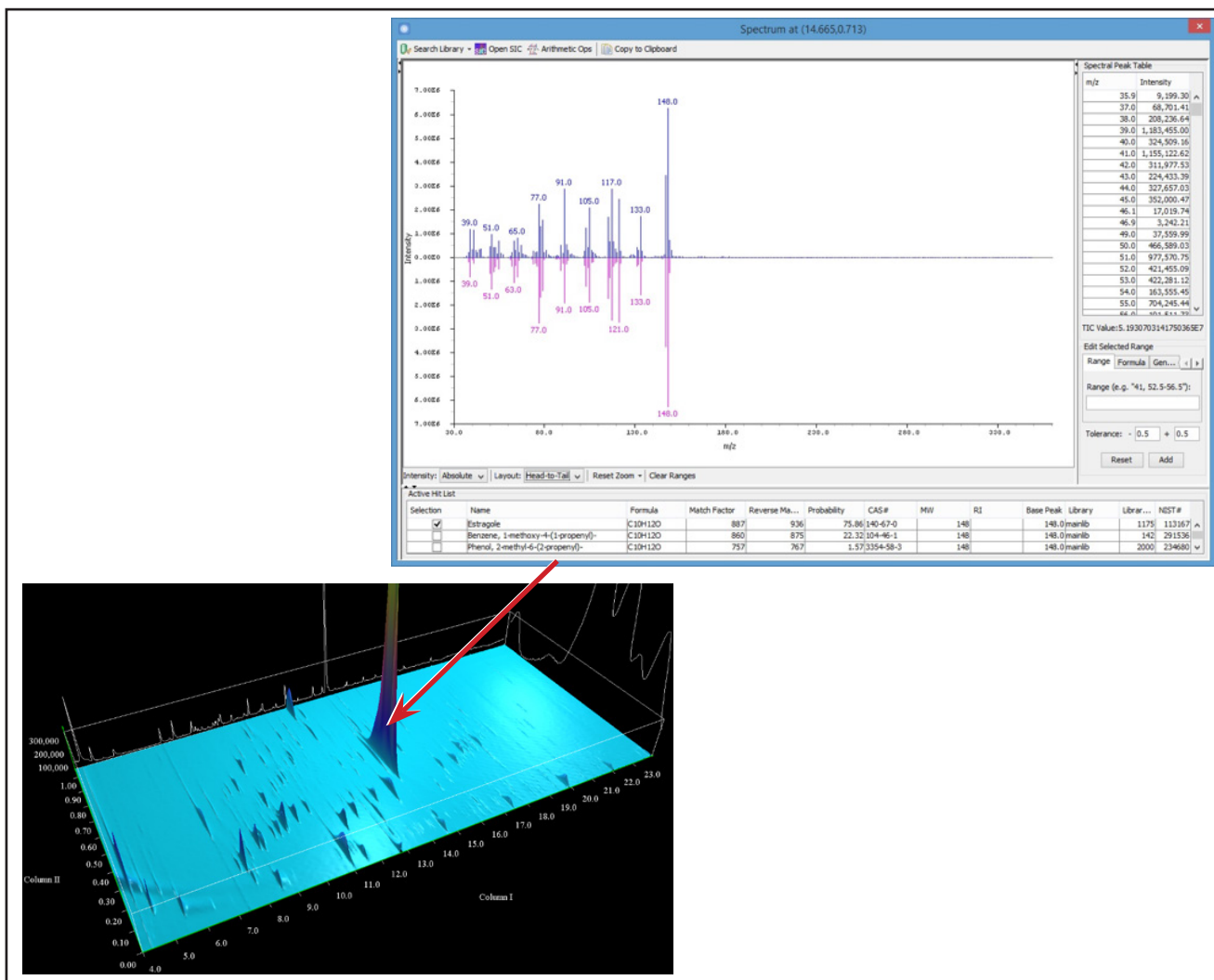


Fig. 5: Estragole peak in sample "Orange Peel, Allspice, Nutmeg"

CONCLUSIONS

- HS-SPME Arrow-GCxGC-MSD/FID is a powerful tool in the chemical profiling of flavours and fragrances
- Reversed flow modulation delivers excellent resolution allowing easier, faster, detailed analysis of samples with complex matrices
- The high sensitivity provided by the Agilent 5977B with HES makes it the ideal tool to compensate for any decrease in sensitivity in GC-MS analysis observed as a consequence of switching from Helium to Hydrogen as carrier gas
- The SPME Arrow technique offers highly sensitive, fully automated, targeted extraction of compounds from complex sample matrices

Legal Statements

CTC Analytics AG reserves the right to make improvements and/or changes to the product(s) described in this document at any time without prior notice.

CTC Analytics AG makes no warranty of any kind pertaining to this product, including but not limited to implied warranties of merchantability and suitability for a particular purpose.

Under no circumstances shall CTC Analytics AG be held liable for any coincidental damage or damages arising as a consequence of or from the use of this document.

© 2019 CTC Analytics AG. All rights reserved. Neither this publication nor any part hereof may be copied, photocopied, reproduced, translated, distributed or reduced to electronic medium or machine readable form without the prior written permission from CTC Analytics AG, except as permitted under copyright laws.

CTC Analytics AG acknowledges all trade names and trademarks used as the property of their respective owners.

PAL is a registered trademark of CTC Analytics AG | Switzerland

Imprint

Date of print: 02.2019

CTC Analytics AG
Industriestrasse 20
CH-4222 Zwingen
Switzerland
T +41 61 765 81 00
Contact: info@ctc.ch

www.palsystem.com

Visit our homepage for more information.